


**A186-2 MASSAGE CHAIR**  
OPERATION INSTRUCTIONS

**SAVE THESE INSTRUCTIONS**

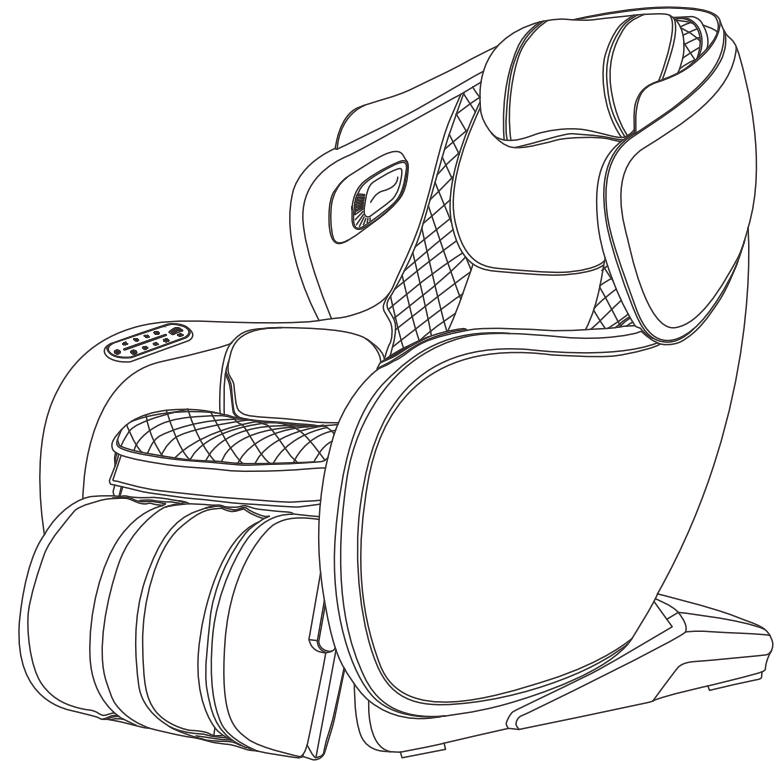

Thanks for purchasing this product. Please read this manual instruction carefully in order to know how to use this product correctly before operation. Please pay more attention to some important safety warning and keep this manual book properly for future use.

Note: Our company reserves the right to change and explain the design of the product. If there is some difference between the product function in market circulation and the description in the manual, it is a normal update of the product, and our company will not inform you otherwise. The specific product function and color shall prevail in kind.

## **CONTENTS**

### **Before Use**

---

Safety Precautions -----	02-04
Place Massage Chair -----	05-06
Before Use -----	07-08
Product Structure -----	09-12

### **Usage Method**

---

Operation Instructions -----	13-17
------------------------------	-------

### **Maintain**

---

Cleaning and Maintenance -----	18
--------------------------------	----



### **Important**

---

Troubleshooting -----	19
Product Specifications -----	20

**Safety Precautions** (must be observed)


- Please read all relevant instructions carefully before use.
- The notes marked here will help you to use the product safely and correctly.
- In order to clearly indicate the extent of harm and damage, precautions are divided into **"Warning"** and **"Caution"** that may be caused by incorrect use. Both of them are important for safety and must be followed.

 <b>Warning</b>	Indicates that serious injury or death may result from improper use.
 <b>Caution</b>	Indicates that misuse may cause minor injury or property damage.


- Please pay attention to the symbol before each safety item.


 <b>Forbidden</b>	Absolute prohibition
----------------------------------------------------------------------------------------------------	----------------------


**The following populations should not use this product**

 <b>Forbidden</b>	<ul style="list-style-type: none"> <li>■ The device is not intended for use by persons (including children) with physical disabilities, sensory or neurological impairments or who lack experience and common sense, unless they are supervised and instructed in the use of the device by those responsible for their safety. Children must be supervised to ensure they do not play with instruments.</li> <li>■ Patients with osteoporosis.</li> <li>■ Patients with heart disease or wearing electronic medical devices such as a pacemaker.</li> <li>■ People who have a fever.</li> <li>■ Pregnant women or women on their menstrual period.</li> <li>■ Patients who are injured or have a skin disease.</li> <li>■ When unsupervised, any children are under the age of 14 and those unconscious are not allowed to use this product.</li> <li>■ People who are told by their physicians to have rest or who feel unwell.</li> <li>■ People whose bodies are wet are forbidden to use this product.</li> <li>■ This product has a hot surface. Persons insensitive to heat must pay attention when using this product.</li> <li>■ If the supply cord is damaged, to avoid danger, it must be replaced by a professional from the manufacturer, its maintenance department or similar department</li> </ul>
-----------------------------------------------------------------------------------------------------	-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------


**Safety Precautions** (must be observed)


<b>Important safety warnings</b>	
 <b>Warnings</b>	<ul style="list-style-type: none"> <li>■ Do not allow children to touch any movable part of this product.</li> <li>■ Use the well-grounded power supply unit suitable for this product.</li> <li>■ Pull out the plug after use or before cleaning to avoid any physical injury or damage to the product.</li> <li>■ Operate this product as instructed in this manual.</li> <li>■ Do not use any parts or accessories which are not recommended.</li> <li>■ Do not use this product outdoors.</li> <li>■ Please read the Operation Instructions carefully before operation.</li> <li>■ Never use this product for any other purpose not listed herein.</li> <li>■ 20 minutes of usage each time is recommended.</li> <li>■ Do not use this product if the cover or leather is broken or damaged.</li> <li>■ Do not drop anything into this product.</li> <li>■ Never fall asleep while using this product.</li> <li>■ Do not use this product if you are drunk or feel unwell.</li> <li>■ Do not use this product within one hour after a meal.</li> <li>■ Do not use this product with excessive force to avoid any injury.</li> </ul>


<b>Safety precautions</b>	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Check whether voltage is consistent with the specification of this product.</li> <li>■ Never insert or pull out the plug with a wet hand.</li> <li>■ Do not let water flow into this product to avoid electric shock or cause damage to this product.</li> <li>■ Do not pull the power cord when plugging or unplugging it. Rough handling is forbidden.</li> <li>■ Do not damage the wire or modify the circuit of this product.</li> <li>■ Do not clean live parts of this product with a wet cloth such as a power switch and plug.</li> <li>■ Be away from this product in case of power failure to avoid any injury if power is restored suddenly.</li> <li>■ Stop using this product at once while it works abnormal, and consult your local supplier.</li> <li>■ Stop using this product if you feel unwell, and consult your health care practitioner.</li> </ul>

<b>Storage conditions</b>	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Storage room temperature: 5 °C to 40 °C.</li> <li>■ Relative humidity (RH) between 20% and 80%; free of corrosive gases, and a well-ventilated room.</li> </ul>


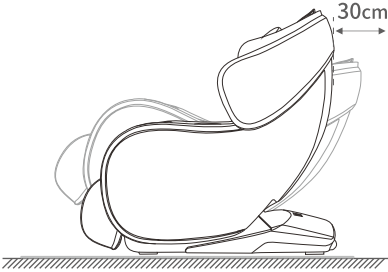
**Safety Precautions** (must be observed)


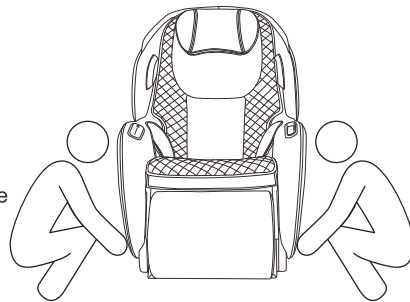
Product service and maintenance	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ This product should be maintained by the service center designated by the manufacturer only. No user should disassemble or maintain it without permission.</li> <li>■ Do not forget to turn off the main power switch after usage.</li> <li>■ Do not use this product if the power socket becomes loose.</li> <li>■ If this product will be left unused for a long period of time, it is necessary to roll up the power cord of this product and store this product in a dry and dust-free environment.</li> <li>■ Do not store this product at a high temperature or near open flame. Avoid prolonged exposure to direct sunshine.</li> <li>■ Please clean this product with dry cloth. Never use thinner, benzene or alcohol.</li> <li>■ The mechanical components of this product are specially designed and manufactured, so no special maintenance is needed.</li> <li>■ Do not use a sharp object to stab at this product.</li> <li>■ Do not roll or pull this product on any uneven ground. Instead, it should be lifted before moving.</li> <li>■ Please use this product intermittently. Do not keep it working continuously for a long period of time.</li> </ul>

Operating environment	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ Do not use this product in a very moist environment such as the bathroom.</li> <li>■ Stop using this product immediately when the ambient temperature changes sharply.</li> <li>■ Do not use this product in any heavily-dusty or caustic environment.</li> <li>■ Do not use this product in a small space or a place without good ventilation.</li> </ul>

Solutions to ordinary malfunctions	
 <b>Caution</b>	<ul style="list-style-type: none"> <li>■ If the motor sounds during operation of this product, this is the normal operation sound.</li> <li>■ If the controller cannot work normally, check whether the plug and the power socket are connected reliably, and whether the power switch is turned on.</li> <li>■ If the rated working time is over, this product will automatically activate the power-off button; If this product keeps operating for a long period of time, the temperature protector will stop this product automatically. The product can be used again after half an hour.</li> </ul>

**Place massage chair**

Installation placement	
<p>Make sure there is enough space so that the massage chair can be tilted, the back end of the massage chair is more than 30cm from the wall.</p> <p>Do not expose the massage chair directly to high temperature environment due to sunlight (such as: in front of heating equipment), and place it to where can avoid high temperature and direct sunlight. It is recommended to lay a soft cushion for long-term use in a designated location. When laying the cushion, the size of the cushion should be enough to cover Keep the range where the machine touches the floor and the range where your feet touch the floor.</p>	 <b>CAUTION</b> 

Product lift instruction	
<p><b>Require team work to lift and move the chair to avoid injury.</b></p> <ol style="list-style-type: none"> <li>1. Please lift the massage chair to move it on vulnerable ground such as wooden floor.</li> <li>2. When moving on a similar wooden floor, it is recommended to lift it. Avoid damage to the casters. When the floor is lifted and moved to the designated position, be careful to keep your hands out of the way to avoid injury, and be careful of your feet being crushed. Do not let go of the unit until the unit is completely flat, the floor may be damaged, therefore, it is recommended to place mats, etc. on the floor.</li> </ol>	 <b>CAUTION</b> 


## Place massage chair

### Move method

Use the rear casters to move the massage chair. Move the chair by tilting it down and backward about 45° as shown in the diagram.

**NOTE:** Excessive force to tilt will cause the chair to tip over completely and product damage.

1. Before moving, please shut down the machine properly to reset the whole machine, and then cut off the power.
2. Please keep the chair vacant when moving. (People, pets or clutter will not be able to stay on the chair.)
3. When moving, do not cross a drop barrier higher than 2mm or a gap greater than 5mm.
4. When using casters to move this product, single continuous movement is not supported; if the distance exceeds 50 meters, please make this product move slowly and evenly.

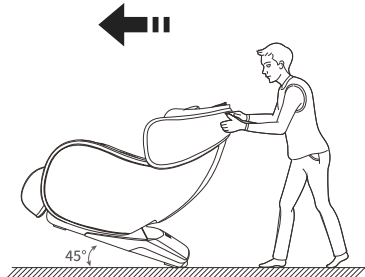


**CAUTION**

## Before use

### Check around


Make sure that there are no people, pets and other objects around the unit.

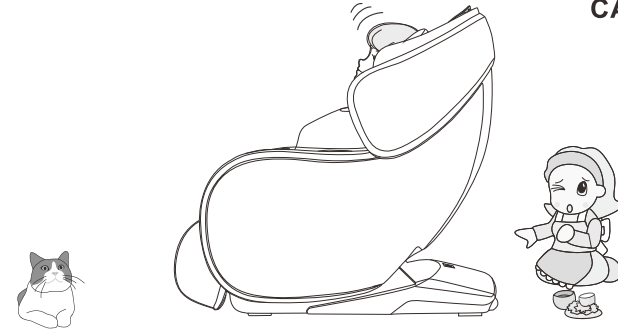

**CAUTION**

### Make sure there are no foreign objects in the gaps of the machine


Before massaging, make sure that no foreign objects are caught in the legs, feet, etc. Under no circumstances allow hands, feet, and head to penetrate into the area between the lower leg mechanism and main body, or the area between the armrest and the capsule. Under no circumstances should hands, feet, or head sink under the fabric or leather covering of the massage area.

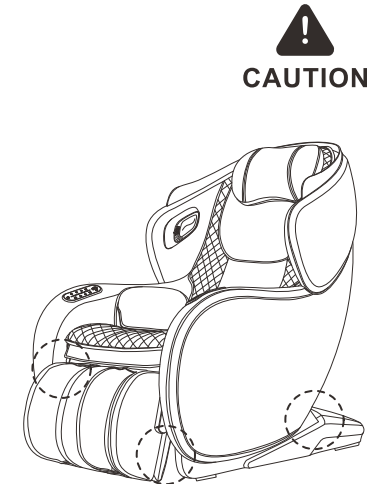


**CAUTION**


**Before use**

### Calf Instructions

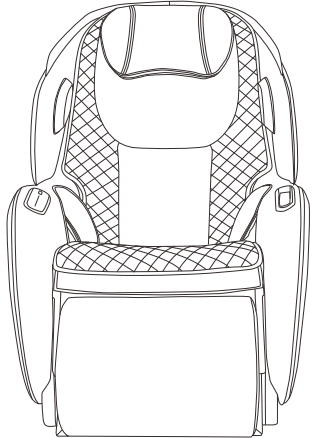

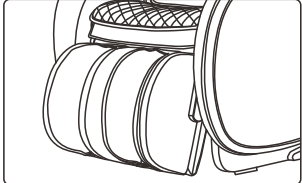

Message mode: flip to enjoy calf air pressure massage

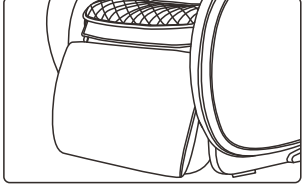


Casual mode: Flip again to become sun stirrups, easy to store, simple and convenient


### Turn on the power

**Caution**


Plug in power cord

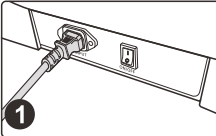

Schematic diagram of power on

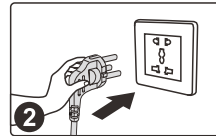

Turn on the power switch ("I" position)

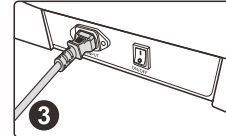


**Warn**

- Please check whether the power cord and power plug are damaged before use.
- Please supply power to this product in accordance with the working power requirements of this product's calibration.
- It is strictly forbidden to use a power supply other than the nominal power supply of this product.
- It is strictly forbidden to connect various power conversion devices to supply power to this product without authorization.


**Product Structure**

**Parts description**

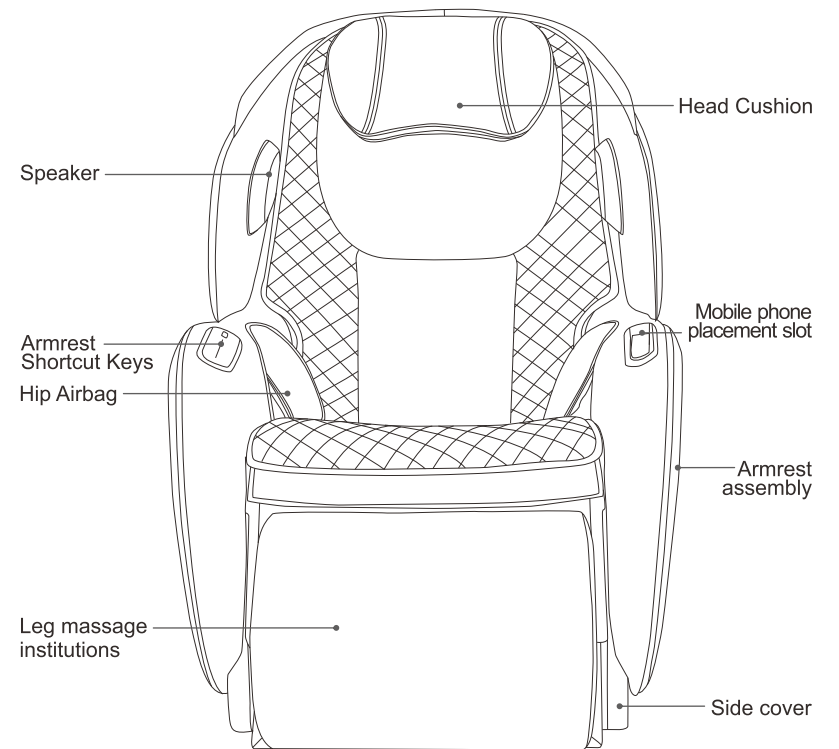


## Product Structure

### Function distribution

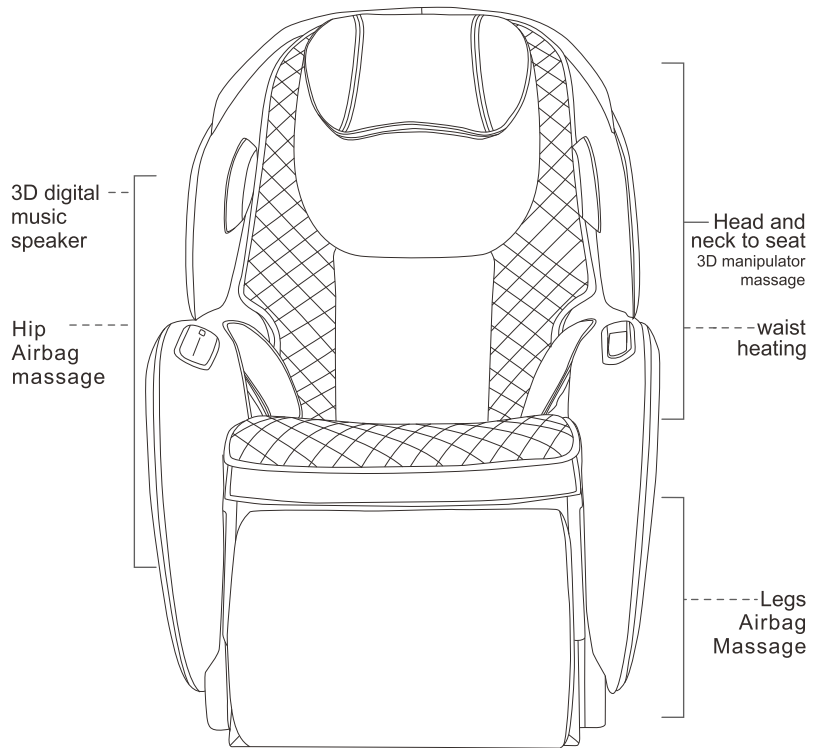





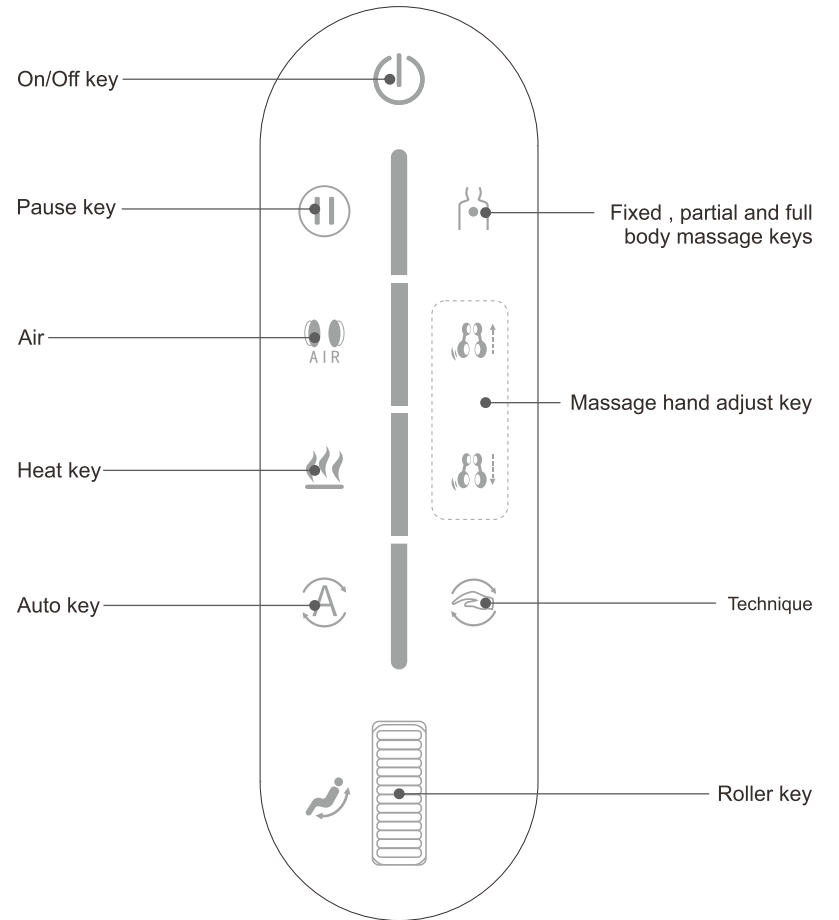





## Product Structure

### I .Right armrest shortcut button instructions





## Product Structure

Function	Description
	Power key: Short press this key to turn ON/OFF the massage.
	Press this key during massage, and all the massage actions will be suspended at once. (i.e. the massage chair will turn off after 20min). Then press this key again to restore all the previous massage actions.
	Air pressure key: Turn on/off air pressure.
	Heat treatment key: Click the key to enable or disable the heat treatment function.
	Auto mode key: each time the key is short pressed, it will be switched to next automatic massage function.
	Fixed, partial and full body massage keys.
	Massage hand adjustment keys: adjust massage acupressure points on the human shoulder. Massage position adjustment (in fixed point and local mode, massage position can be adjusted).
	
	Technique: For details, please refer to page 14 of this manual.
	Roller key: control backrest up or down, see page 15 of the manual for details.









## Operation Instructions

### II .Start massage

- ①. Press the power button to start the massage, the default is relax massage.
- ②. Automatic angle adjustment (you can press the "A" button to skip).
- ③. Automatic human shoulder position detection, keep the sitting posture fit with the massage chair (you can press the "A" button to skip).
- ④. Shoulder position fine-tuning, if the detected position deviates from the actual shoulder position, adjust it to the appropriate position through the up and down buttons of the shoulder position (you can press the "A" button to skip), there are 11 shoulder positions in total. If there is no adjustment within 10 seconds, it will enter into the massage function.

Function	Description
	Move the massage balls upward to a suitable position
	Move the massage balls downward to a suitable position

#### • Automatic massage












Function	Button	Description
Automatic massage	 White	Full body massage: deep massage to rejuvenate the body.
	 Orange	Comfort: The most comfortable massage experience.
	 Red	Shoulder & neck: focus on shoulder and neck massage.
	 Blue	Back&waist: focus on back and waist massage.
	 Mauve	Relax: the depth of the soft and comfortable massage, relieve body fatigue.
	 Dark purple	Spinal massage: Focus on spine massage and fully relieve the fatigue and pain from your body.
	 Green	Beautiful buttocks shaping: effective massage for the buttocks and thighs, plays a role in shaping the body, reducing the accumulation of excess body fat.
	 Yellow	Leisure Heath: for sedentary sitting in the office, massage can increase the amount of physical activity and improve the body's function.

Note: Automatic program 1, 2, 3 contains background music, long press "A" for 3 seconds to turn on/off the background music.


When the power is on, press and hold "☺" for 3 seconds to turn on/off the LED light.

## Operation Instructions


### • Manual massage

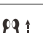Function	Button	Massage ball adjustment	
Technique	 White	Knead.	
	 Light blue	Beat.	
	 Red	Knead pat.	
	 Green	Shiatsu.	
	 Orange	Knock.	
	 No	No action.	
Back stretch	 Red	Fixed position and spot massage.	
	 Purple	Small area local round trip massage.	
	 White	A round trip massage.	
		Message position adjustment (in fixed point and local mode, massage position can be adjusted).	Manipulator uplink
			Manipulator descending


### • Pneumatic massage

Function	Button	Description
pressure		Press this key to turn the air pressure on or off.

### • Heated massage


Function	Button	Description
Heat		Press this key to turn on or off the back far infrared heat therapy, open 3 minutes after the warmth.

### • Pauselstart massage button

Function	Button	Description
Pause/start		Press this button in the process of massage, and all massage movements will be suspended urgently when time continues to countdown), and then press this key to resume all before suspension massage.

## Operation Instructions

### • Adjust lying Angle

Function	Button	Description
Sit up		Turn up the roller, the backrest frame slowly rises, the small leg frame slowly declines, and emits "beep... Beep..." Release this key to stop up immediately.
Lay down key		Turn down the roller, the backrest frame slowly declines, the small leg frame slowly rises, and emits "beep... Beep..." Release this key to stop lying down immediately.

### III. Operation Instructions for Bluetooth Speaker Play

1. Enter the smart device main menu, Open "Settings" and select the "Bluetooth" options. Open "Bluetooth", and Bluetooth will automatically "Search for a New Bluetooth Device".
2. The smart device will search for the corresponding Bluetooth model (e.g. IMCM-XXXX), Click the "Match" button, if you are asked whether to match with the device, you just select "Match" to complete matching.
3. Turn on "Music Player" in the smart device and select the song to play.

### IV. Voice Control Mode

- When the massage chair is turned on or in standby mode the sound "Hi Alice" or "Hey Alice" will be produced to enable the voice control function. Alternatively, you can press the voice control button on the armrest to enable the voice control function.
- After such function is enabled successfully, the massage chair will reply "I am here" At this moment, once a corresponding voice control entry to be implemented is given within six seconds, the corresponding voice control entry to be implemented is given within six seconds, the corresponding voice control function can be performed after a reply is received from the massage chair if no correct voice control entry is given within six seconds, the voice control function automatically will be off automatically and needs to be enabled again.

## Operation Instructions


· Note: When using the voice control function, the speech speed should not be too fast, and the voice commands should be issued according to the voice control entries. After each successful use of the voice-activated entry, it is recommended to give out the voice command after 1-2 seconds.

Voice Control Entry		
Entry: "Hi Alice" or "Hey Alice". Answer: I am here.		
No.	Voice control entry	Answer entry
1	Massage on	Answer: Ok, massage on
2	Massage close	Answer: Ok, massage close
3	Comfort Massage	Answer: Ok, comfort massage
4	Relax massage	Answer: Ok, relax massage
5	Full body massage	Answer: Ok, full body massage
6	Neck and shoulder massage	Answer: Ok, neck and shoulder massage
7	Back and waist massage	Answer: Ok, back and waist massage
8	Open the air pressure	Answer: Ok, open the air pressure
9	Close the air pressure	Answer: Ok, close the air pressure
10	Up the seat position	Answer: Ok, up the seat position
11	Down the seat position	Answer: Ok, down the seat position
12	Change the other mode	Answer: Ok, change the other mode
13	Go little down	Answer: Ok, go little down
14	Go little up	Answer: Ok, go little up
15	Fixed massage	Answer: Ok, Fixed massage
16	Stop adjustment	Answer: Ok, Stop adjustment


## Operation Instructions

### V. End the massage


1. During the massage process, press the power key off touch screen controller to immediately turn off all massage functions, the backrest and leg rest will reset, or when the massage timer expires, it will automatically turn off all massage functions, but the backrest and leg rest will not reset.

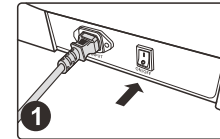

Power off ("0" position)

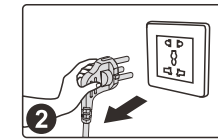

Take off the power plug

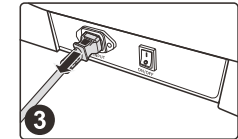


Take off the product power cable.

## Cleaning and Maintenance

### Clean PU leather and plastic parts

Wipe the leather with a soft, dry cloth.

- Do not use chemicals such as thinner, gasoline, alcohol, etc. to clean the unit.
- When using commercially available leather care products (rags), please follow the instructions for use.
- If the PU leather is particularly dirty, wipe it by the following methods.
  - ① Immerse a soft cloth in water or hot water containing 3% to 5% neutral detergent, then wring it out thoroughly.
  - ② Use the above soft cloth to wipe the leather surface.
  - ③ Wipe off the cleaning solution with a cloth soaked in clean water and thoroughly wrung out.
  - ④ Wipe with a wrung soft cloth.
  - ⑤ Let it air dry naturally.
- When it is difficult to remove stains, please immerse commercially available melamine foam in neutral detergent, and then wipe the machine with it.
- Do not use a hair dryer to dry the surface.
- As the color of the fabric will stain the surface of the PU leather, be careful when wiping the chair with denim and colored fabrics.
- Do not keep these parts in contact with plastic for a long time, as this will cause discoloration.
- PU leather may be discolored, so when using hair dye, please cover the touched part with a towel.

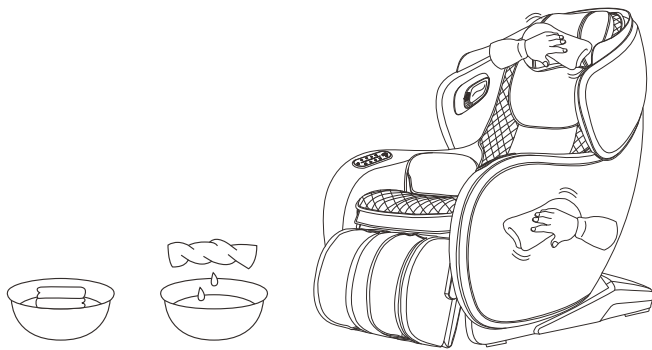

## Troubleshooting

<b>Fault</b>	When the calf mechanism or the backrest reaches a certain position, press the corresponding adjustment button and it will no longer change, and a continuous "beep, beep" sound will be emitted.
<b>Solution</b>	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.
<b>Fault</b>	The function does not work or some functions do not work after the machine is turned on.
<b>Solution</b>	After the machine is turned on by pressing the switch key, it is necessary to manually select an automatic massage program or other working state. If it is not operated, the machine will automatically shut down after 20 minutes. If the manual operation mode is selected, the functions that are not selected may be disabled, and the working mode needs to be manually set in the manual mode.
<b>Fault</b>	The machine makes a leather rubbing sound when adjusting its posture.
<b>Solution</b>	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.
<b>Fault</b>	The massage wheel did not reach the shoulders or neck.
<b>Solution</b>	If the head does not touch the head cushion or the back does not touch the backrest, the shoulder position may be lower than the actual position during the body shape detection process. Sit in the deepest position of the seat, rest your head on the headrest, and start over from the beginning.
<b>Fault</b>	The unit is damaged. The power cord or power plug is abnormally hot.
<b>Solution</b>	To prevent accidents, be sure to contact an authorized service center.

Please feel free to contact the designated service personnel for any questions or doubts about this product.

## Product Specifications

---

Model: A186

Configure: 2

Name: Massage Chair

Rated voltage: 110-120V~ 60Hz

220-240V~ 50Hz/60Hz

Rated power: 80W

Rated time: 20min

Safety structure: Class I

Material description: PVC, PA, steel parts, electric and electron components

MATERIAL	PU	LEATHER	CLOTH	WOOD